

Fotona[®]
choose perfection

60
1964
2024

years of innovation



One Device, Multiple Solutions

StarFormer[®]

Magnetic Body Sculpting and Muscle Strengthening

STARFORMER®



StarFormer® is Fotona's latest innovation in body shaping. This breakthrough technology utilizes magnetic pulse energy to stimulate muscle contraction. As a result, muscles strengthen, tone and firm.

The system provides two treatment options:

- **TightWave®**
muscle strengthening, toning and firming of weakened areas
- **IntimaWave®**
strengthening of pelvic floor muscles to treat stress urinary incontinence

Unique Quattro modality


Enabling simultaneous action of 4 applicators and maximizing your clinic's productivity.

How it Works

StarFormer® uses magnetic stimulation technology to selectively target muscle tissues. Magnetic pulses harmlessly penetrate into the muscle layers causing them to automatically contract. The end result is increased muscle mass.

Extra Large Applicator

- 23.5 cm diameter
- The preferred choice for gluts and abdomen

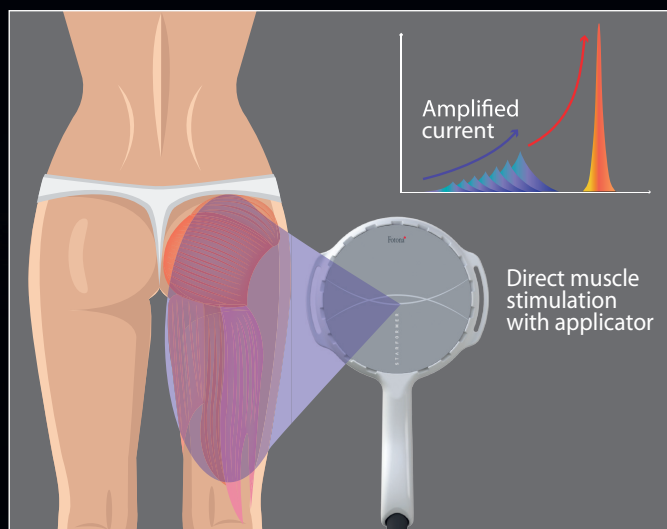

Magnetic body sculpting and muscle strengthening is a breakthrough treatment for effectively strengthening, shaping and firming different areas of the body:



Convenient Applicator Holder

- For easier application e.g. when treating biceps and triceps
- Alternative to straps

- Glutes
- Abdomen
- Hamstring area
- Arms
- Back muscles

The treatment is performed with the patient lying down fully clothed. One session generally lasts between 20 to 30 minutes with patients resuming their daily activities immediately after treatment.


TIGHTWAVE®

INTIMAWAVE®

StarFormer® can additionally be used to stimulate pelvic floor and lower back muscles, utilizing two applicators to treat all types of male or female urinary and faecal incontinence. During treatment the patient remains fully clothed and sits on a comfortable IntimaWave® chair. The magnetic applicators in the chair stimulate muscle contraction in the entire pelvic floor area, increasing the strength and endurance of the pelvic floor. The result is effective and long-term relief of uncontrolled urine flow.




Before

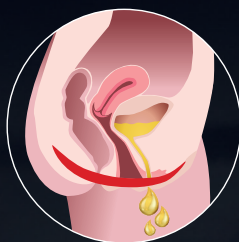

Pelvic floor muscles insufficiently support pelvic organs, affecting bladder control.

After

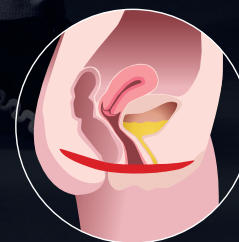


Magnetic stimulation induces pelvic floor muscle strengthening and the return of normal continence function.

STARFORMER®

Key Benefits

- Safe, quick and easy muscle strengthening
- Non-invasive
- Improvement after only a few sessions
- Walk-in/walk-out procedure
- High success rate and patient satisfaction

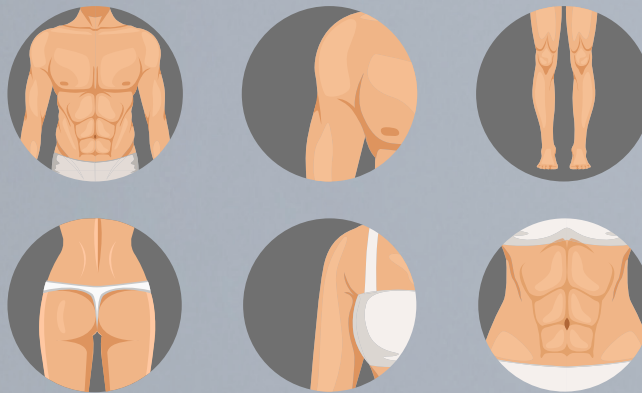

Device name:	Fotona StarFormer®
Input voltage:	100-240V AC, 50~60Hz
Power:	1500 VA
IP protection:	IPX0
Size (W x L x H):	421 mm x 630 mm x 843 mm
Weight:	50 kg
IEC 60601/1 classification:	Class I, Type BF
Regulation (EU) 2017/745 classification:	Class IIa
Available accessories:	IntimaWave Chair
	Applicator Radial Large
	Applicator Radial XLarge

CE
2265

All Fotona medical devices are CE marked and approved to be sold in the EU. For countries where specific national approvals or clearances are required, some of the products and/or applications may not yet have been approved. Please check with Fotona, your local Fotona distributor or your national regulatory body about whether a specific product or application has been approved to be marketed and sold in your country.

For related patents see: www.fotona.com/patents

www.fotona.com

Fotona, d. o. o.
Stegne 7
1000 Ljubljana
Slovenia
EU

Fotona, LLC
4343 W Royal Lane #116
Irving, TX 75063
USA

Fotona Beauty Light, (Suzhou)
Medical Devices Co, Ltd.
No 2, Zengfu Road, Guli Town
Changshu City, Jiangsu Province
CHINA, 215515

Fotona France SARL
47 Boulevard de Courcelles
75008 Paris
France
EU

Fotona GmbH
Hohlbachweg 2
73344 Gruibingen
Germany
EU

Fotona Japan, Ltd.
Tokyo
Japan